

◆◆◆ **Dr.Shine®**

by  **healthshine**
LIVING IN GOOD HEALTH



AirVista BPAP-30 PRO

For Sleep Disorder Patients

Air Vista BPAP Pro is design to deliver advance non -invasive ventilation with precision , comfort, and reliability. For patients requiring higher pressure support, it combines cutting edge technology with user friendly features to ensure effective therapy and restful nights.



EU Certified

Real-Time Data Monitoring

Track the monitoring treatment data
in real time

3.5" Colorful LCD screen

Provides vivid information in one page and
more legibility; User interface design to
enable fast-setup.

Master Core-Technology

Imported Swiss pressure sensor
cooperates with innovative
algorithm
to fulfill greater sensibility.



Dial Design

Rotate one key to
complete
various operations

Lower dBA


Less noise achieved by
high-quality blower

Heating Humidifier

5 temperature levels for humidifier are available;
Designed to prevent back flow and detachable
to clean easily.



Leak Compensation

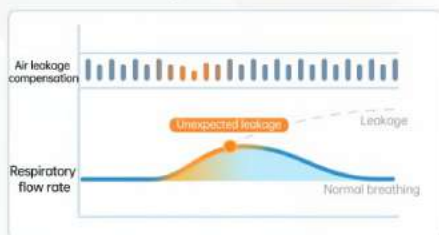

Main Interface



Multiple Respiratory Modes

CPAP, S, T, S/T, APCV suitable for various respiratory insufficient patients

Auto-Sensibility

Isens and Esens are automatically adopted to individuals treatment

VAF

Maintains a target tidal to meet patient's needs

Variety of Alerts

Safer use is ensured by alerts, including power down, air leakage, high/low pressure, etc.

InspTime Limitation

User can choose a period of time to inhale by setting the max- and min-InspTime, which ensures effective treatments.

Technical Specification	
Mode	Cpap, S/T, S, T, APCV
Display Screen	3.5" Color LCD
Ipap Pressure	4-30 Cm H ₂ O
Epap Pressure	4-25 Cm H ₂ O
Isens	Auto, 1-6 Levels
Esens	Auto, 1-6 Levels
I Slope	1-6 levels
I Slope Time	0.5 - 4 sec
Max. Insp. Time	0.5 - 4 sec
Min. Insp. Time	0.5 - 4 sec
BPM	3-40 bpm
Target VT	50-2000ml
Ramp Time	0-60min
VAF Setting	Yes
Auto On	Yes
Auto Off	Yes
Altitude Compensation	Up to 8,200 Feet
Leak Compensation	Up to 150 L
Heated Humidifier	0-5 levels, 50-70 °C
Noise Level	≤30dB
Configuration	
Data Storage	365 days, Micro SD card
PF Curve	Yes
Dimension	280*200*112 mm
Weight	2.2 kg
Power Supply	AC 220 V, 50Hz
Power Consumption	≤80 Watt with Humidifier
Alarms	Power Failure, High & Low Pressure, Low MV, Low Voltage, High Leakage, Airway Block, Humidifier Failure, Memory card Full
Warranty	2 years on device