

VK300

The Vital Key of Life with Intelligence and Innovation



Better Vision, Better Performance

- 15.6" touch screen with navigation knob for easy operation
- The big font view and simultaneous display make it easy to observe
- Intuitive user interface ensuring efficient use
- Ergonomic designed trolley compatible with air compressor

Quick
Start

High
Performance

Easy
Maintenance

Invasive to
Non-Invasive

Pediatric
to Adult




More Value, More Possibility


Minimum tidal volume at 10ml

12+ ventilation modes

Advanced Ventilation Modes: BiPPV, APRV, PRVC, etc.

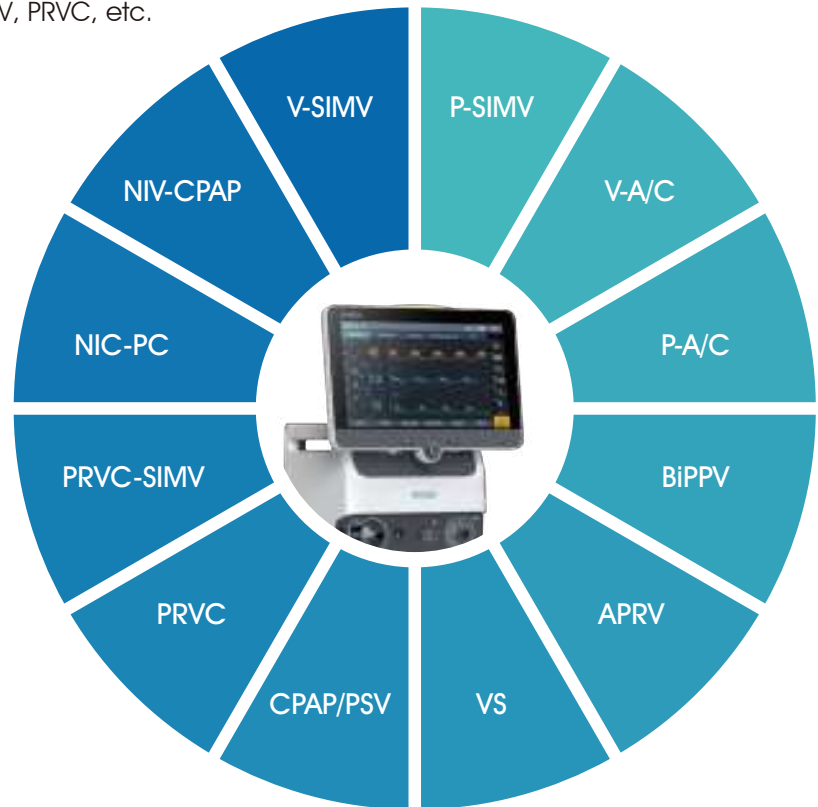


Access to Easy Airway Management

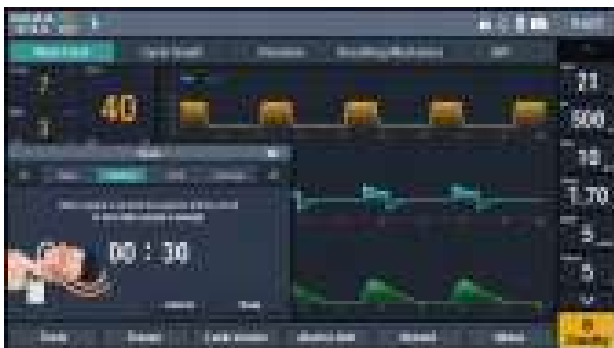

Suction Maneuver

- Reduce suffocation by preoxygenation
- Coordinate with suctioning
- Ensure no interruption during the procedure

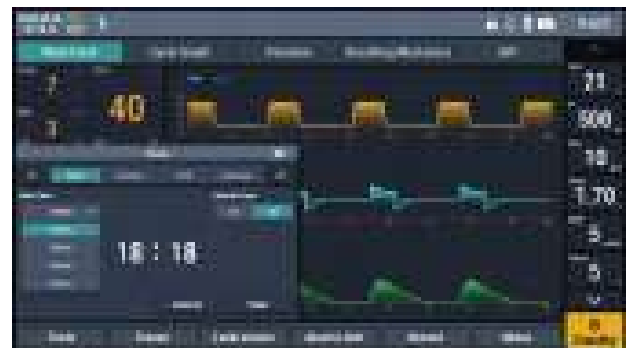


Medication Nebulization

- More easy absorption of the drug
- Select-able neutralizing time
- Improve airway tensity effectively




Practical Clinical Tools

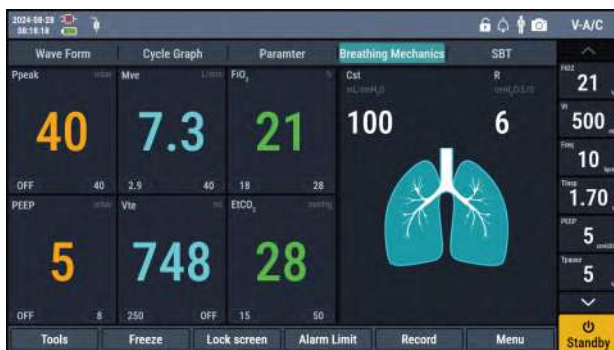

Dynamic Lung

Vivid graphic display
Stress Index
Dynamic and static compliance

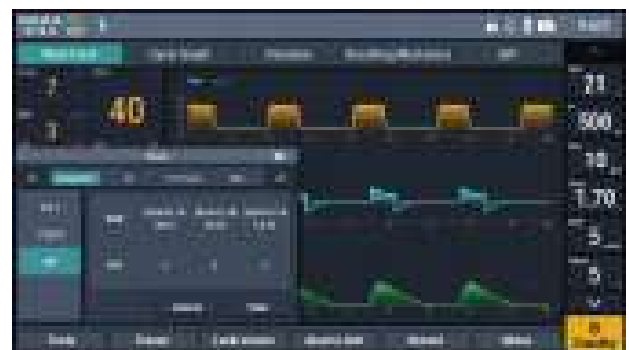

Weaning Parameters

Weaning parameters include P0.1,
Negative Inspiratory Force (NIF),
Rapid Shallow Breath Index (RSBI), etc.

Lung Recruitment



One-key
activation



Better
oxygenation

P-V Tools




Automatic
measurement



Automatic
analysis



Easy Management, Escort Every Life

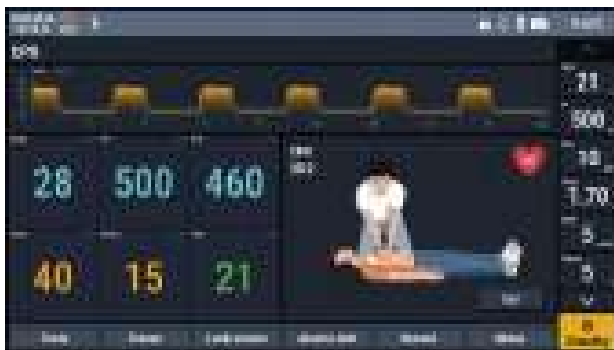


CPR

- Automatic adjustment of pressure limit
- Avoid false triggers and confrontation
- Voice guidance for more effective rescue

Compensation

- Lung compliance compensation
- Automatic tube resistance compensation
- Available for invasive & non-invasive ventilation

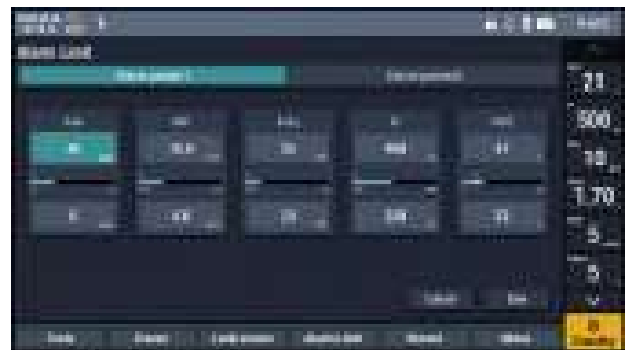

Alarms System

- 360° alarm lamp
- Visualized and audible alarm
- Traceable events

EtCO₂ Module (optional)

- Real-time CO₂ measurements
- Better user experience using a non-invasive method

Global Markets



Amoul[®] Ambulanc (Shenzhen) Tech. Co., Ltd.



📍 Building #8, Skyworth Innovation Valley Phase II, TangTou Road No.1,
ShiYan, Bao'An District, ShenZhen, China
☎ +86 755-26072286
🌐 www.amoulmed.com
@ info@amoulmed.com

